




An ISO 9001 & 14001 Co.

## HOSPITAL & MEDICAL EQUIPMENT

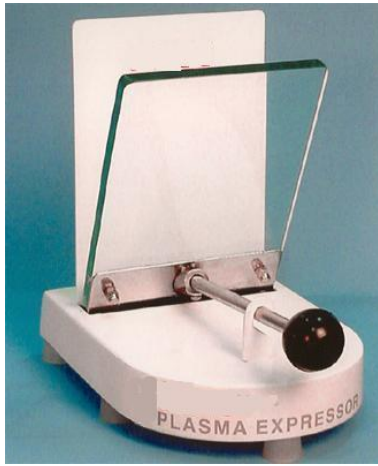

### Plasma Extractor

Item Code: BB-05

Rajas **Plasma Extractor** is designed to extract blood component from centrifuged bags. The spring loaded front panel applies pressure on the collected bag causing the liquid to contain in transfer bag. Its sturdy construction makes it durable and reliable, easy to use and portable. The unit is available in both versions automatic and manual. Sensor infra red, clamping control motor activated, pressure plate acrylic, alarm audio visual

## RAJAS ENTERPRISES (INDIA)

Nanhera Chowk, Nanhera, Post Office Kuldeep Nagar  
Ambala Cantt- 133004 (Haryana) India  
E-mail : [rajasmedical@gmail.com](mailto:rajasmedical@gmail.com)  
[www.rajasindia.com](http://www.rajasindia.com)